

MIRACLE SMILE

WATER FLOSSER

DELUXE PRO

FLOSS IN SECONDS!+
FAST • EASY • EFFECTIVE

**INSTRUCTION
MANUAL**




Thank you for purchasing the Miracle Smile® Deluxe Pro water flosser, the easy way to floss in seconds!

Customer Satisfaction is our #1 Priority,
but if you are not satisfied for any reason, please contact us directly.


MiracleSmile@rephelpdesk.com

WHAT'S INCLUDED

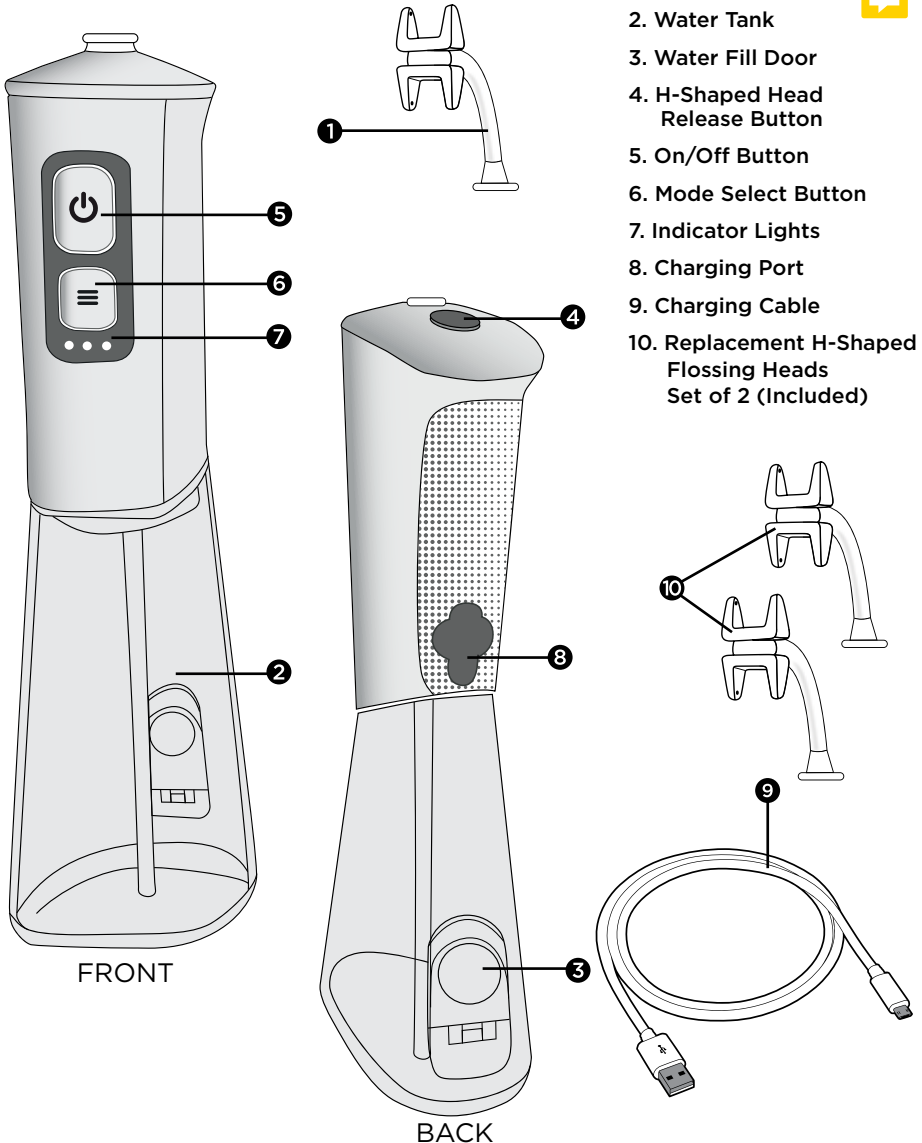


HOW TO CHARGE

We recommend you fully charge your Miracle Smile® Deluxe Pro before the first use.

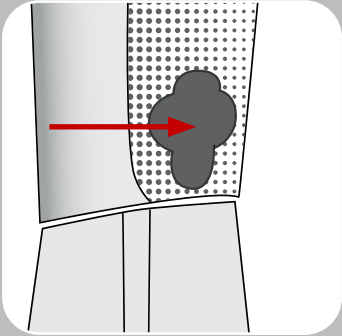

1. Open the protective silicone cover on the Charging Port on the back of the unit.

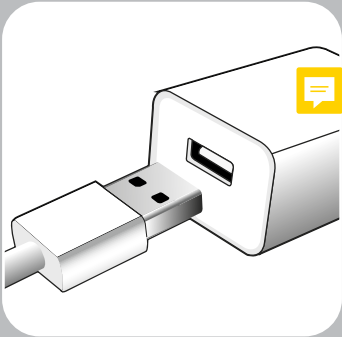


2. Insert small end of Charging Cable into Charging Port and large end into a USB port or USB power adapter (not included).




3. Indicator Lights will blink when charging. Unit is fully charged when all three lights stop blinking and remain fully lit. Remove the charging cable from unit.


HOW TO USE

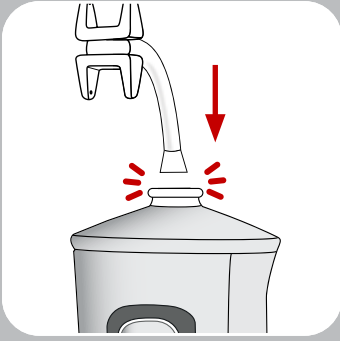

1. Insert H-Shaped Flossing Head into top of unit and press until it snaps into place.

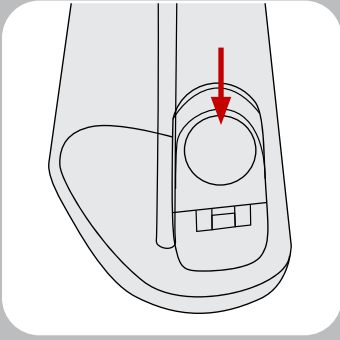


2. Open the Water Fill Door on the back of the Water Tank.

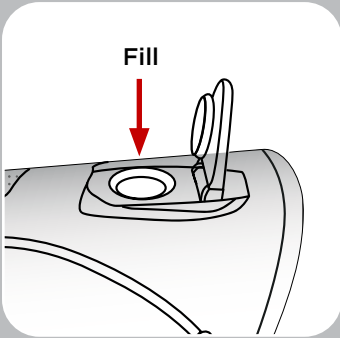


3. Carefully fill Water Tank with water until full (7 ounces).

NOTE: For an enhanced experience, add your desired amount of mouthwash to Water Tank before filling it completely.

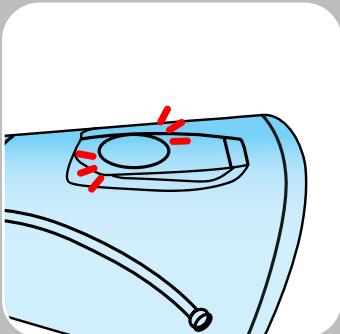


4. Close Water Fill Door until it snaps into place securely to avoid leakage.

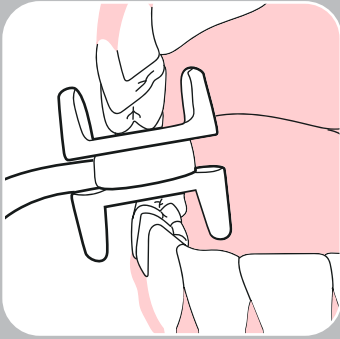


5. Insert H-Shaped Flossing Head into your mouth. Place bottom of H-Shaped Flossing Head around bottom teeth and slowly close your mouth to position top teeth into top of H-Shaped Flossing Head.

NOTE: Press ON/OFF Button after inserting into your mouth.


Select between Low, Medium and High pressure settings with Mode Select Button.

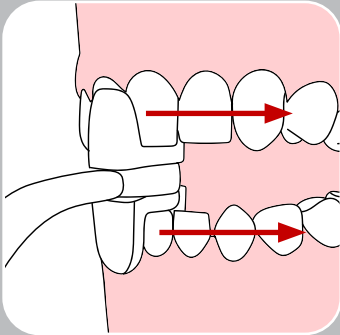


6. Gently glide H-Shaped Flossing Head across teeth from  side to the other in a contiguous, controlled motion. If needed, refill Water Tank and repeat for your desired results.




7. When done, turn off by pressing On/Off Button.

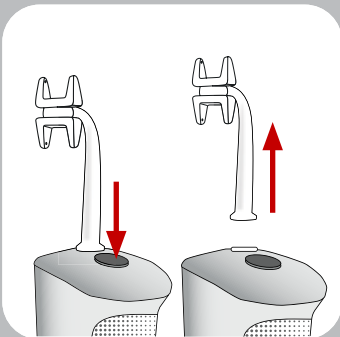

8. To remove H-Shaped Flossing Head, press H- Shaped Flossing Head Release Button and gently pull off unit.

MAINTENANCE

After each use,  empty water tank, clean any leftover debris, and let the unit air dry. Store in a cool, dry place until next use.

To remove/replace H-Shaped Flossing Head, press the release button on the top of unit, and pull out until removed.

WARNINGS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following: **READ ALL INSTRUCTIONS BEFORE USING.**

WARNING: To reduce the risk of injury to persons:

- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not direct water into the nose or ear.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use tips and accessories that are included.
- Do not drop or insert any foreign object into any opening on the product.
- Do not use this product outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce performance and will shorten the life of the product.
- Fill reservoir with warm water or other dental professional recommended solutions only.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.

WARNINGS (CONT.)

- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- Close supervision is necessary when this product is used by, on or near children or invalids.
- Children should be supervised to ensure that they do not play with the appliance.



©2023 Ontel Products Corp
Distributed by Ontel Products Corp.,
Fairfield, NJ 07004,
under license from IMG Int'l, Ltd .
Scarsdale, NY 10583 and Spa Dent Inc.,
Waterloo Canada N2K3T6

Miracle Smile® is a registered trademark of Spa Dent, Inc.,
Waterloo Canada N2K3T6,
used under license.

Made In China
Patents Pending
MCSML-IM-01-0723-300



Due to continued product enhancements,
product image may differ slightly from actual product purchased.